

SorbSan® SILVER Alginate Wound Dressing



PRODUCT DESCRIPTION


SorbSan® SILVER Alginate Wound Dressing is an antimicrobial wound dressing for the management of moderately to highly exuding chronic and acute wounds where there is an infection or an increased risk of infection.

SorbSan® SILVER Alginate Wound Dressing is a highly absorbent, sterile, single use, nonwoven pad or ribbon composed of calcium alginate and 0.95% ionic silver. The dressing absorbs high amounts of wound fluid and the microorganisms contained in the wound fluid. The absorbed wound fluid creates a soft and cohesive gel which assists in maintaining a moist environment and aids in the removal of non-viable tissue from the wound (autolytic debridement). The moist wound healing environment and control of dressing microorganisms supports the body's healing process and may help reduce the risk of wound infection.



DETAILS

SorbSan® SILVER Alginate Wound Dressing is sterilized by Gama Radiation (R). SorbSan® SILVER Alginate Wound Dressing is CE marked under EU Directive 94/42/EEC, and is available for distribution globally.

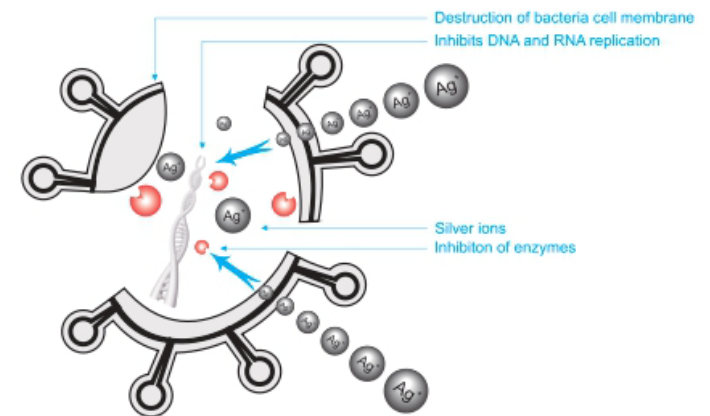

CHARACTERISTICS

- 1 Antimicrobial efficiency for 7 days
- 2 High absorption and fluid retention
- 3 Wet integral and easy application and removal
- 4 Reduced pain during wear and at dressing removal
- 5 Effective against a wide spectrum of micro-organisms

SorbSan® SILVER Alginate Wound Dressing

INDICATIONS

SorbSan® SILVER Alginate Wound Dressing can be used in the management of moderate to heavily exuding chronic wounds including diabetic wounds, venous wounds, pressure ulcers, cavity wounds, etc.

SorbSan® SILVER Alginate Wound Dressing is packed in foil pouches, then into cartons and shipper cases. The packaging configuration is shown in the table below:

Regular Sizes	Qty per box	Qty per shipper case
5x5cm	10pcs/box	10 boxes/shipper case
10x10cm	10pcs/box	10 boxes/shipper case
10x20cm	5pcs/box	10 boxes/shipper case
2.5x45cm	5pcs/box	10 boxes/shipper case
More sizes is also provided on request		



Foshan United Medical Technologies Ltd

We combine our expertise in advanced wound and skin care products, and our manufacturing and production knowledge, to produce superior performing products for our global B2B customers.

Tel: +86-139-28979836
Fax: +86-757-81207311
Email: sales@baihemedical.com
Website: www.baihemedical.com

UMT